

For Research Use Only

# RAI1 POLYCLONAL ANTIBODY



Catalog Number: 17084-1-AP

## Basic Information

Catalog Number:  
17084-1-AP

Size:  
31 µg/150 µl

Source:  
Rabbit

Isotype:  
IgG

Purification Method:  
Antigen affinity purification

Immunogen Catalog Number:  
AC9046

GenBank Accession Number:  
BC021209

GeneID (NCBI):  
10743

Full Name:  
retinoic acid induced 1

Calculated MW:  
1906aa, 203 kDa; 966aa, 104 kDa

Observed MW:  
250-260kd

Recommended Dilutions:  
WB 1:200-1:2000

## Applications

Tested Applications:  
ELISA, WB

Species Specificity:  
human, mouse, rat

Positive Controls:

WB : HeLa cells;

## Background Information

### Notable Publications

Author	Pubmed ID	Journal	Application
--------	-----------	---------	-------------

## Storage

Storage:

Store at -20°C. Stable for one year after shipment.

Storage Buffer:

PBS with 0.02% sodium azide and 50% glycerol pH 7.3.

Aliquoting is unnecessary for -20°C storage


For technical support and original validation data for this product please contact:

T: 1 (888) 4PTGLAB (1-888-478-4522) (toll free in USA), or 1(312) 455-8498 (outside USA)

E: [proteintech@ptglab.com](mailto:proteintech@ptglab.com)  
W: [ptglab.com](http://ptglab.com)

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

## Selected Validation Data

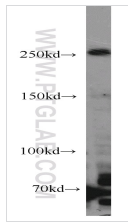

HeLa cells were subjected to SDS PAGE followed by western blot with 17084-1-AP(RAI1 antibody) at dilution of 1:500