

NR1D2 Polyclonal Antibody

Catalog number: 13906-1-AP

Size: 30 µg/150 µl


Source: Rabbit

Isotype: IgG

Synonyms:

NR1D2; BD73, EAR 1R, NR1D2,

Rev erb beta, Reverb Beta, RVR

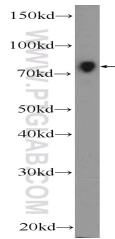

HeLa cells were subjected to SDS PAGE followed by western blot with 13906-1-AP(NR1D2 Antibody) at dilution of 1:300

Applications

Tested applications:	ELISA, WB
Species specificity:	Human,Mouse,Rat; other species not tested.
Calculated NR1D2 MW:	579aa,65 kDa
Observed NR1D2 MW:	70 kDa
Positive WB detected in	HeLa cells, HEK-293
Recommended dilution:	WB: 1:200-1:1000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag4881
GenBank accession number:	BC045613
Gene ID (NCBI):	9975
Full name:	Nuclear receptor subfamily 1, group D, member 2

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.