Cole-Parmer®

Cole-Parmer® WB-200 Series Water Baths

- Digital display with timer and temperature settings
- Corrosion-resistant 304 stainless steel tank
- Built-in overtemperature protection
- Supplied with stainless steel gabled cover
- 2-year warranty



Cole-Parmer® WB-200 Series Water Baths

Use these general-purpose digital water baths for a wide range of laboratory applications. Easily set and view current and set temperatures, timer, and bath status on the large LCD. In the event of a fault, the overtemperature protection system activates the alarm mode turning off all heating elements and lighting a red alarm indicator on the screen. Bath automatically resumes operation once the temperature drops back below the overshoot threshold.

Stainless steel interior and powder-coated steel exterior offer excellent durability. Attached 304 stainless steel gable cover reduces heat loss and ensures uniform temperatures. The built-in timer has a range of 0 to 5999 minutes; delay heating feature can be set to start the cycle when it is convenient for you. An audible alarm sounds intermittently to indicate end of heating cycle.

Perforated bottom plate (included) prevents direct contact with heating elements. Use the optional tube racks to hold various tube sizes in the bath. The number of module spaces listed (see table) determines how many trays fit in each bath. Knobbed handle makes it easy to place or remove the rack. Use the optional holed plates to accommodate flasks in the bath; simply remove the rings to fit flasks from 33 to 92 mm in diameter.

Order optional tube racks and holed plates separately (see next page).





Perforated bottom plate (included)

| Capacity | Temperature range | Temperature stability | No. of module widths accepted | Bath opening (W x H x D) | Overall dimensions (W x H x D) | Item weight | Heater | Electrical (VAC, Hz) | Item number |
|---------------------------|-------------------------|--------------------------|-------------------------------------|--|---|------------------|--------|-------------------------|----------------|
| 6.6 L | Ambient +5 | ±0.2 °C | 2* | 11 ¹³ / ₁₆ x 6½ x 5 ¹⁵ / ₁₆ in | 18 ¹⁵ / ₁₆ x 8 ¹¹ / ₁₆ x 13 ³ / ₈ in | 20 lb (9 kg) | 600 W | 120, 60 | 78905-20 |
| 0.0 L | to 100 °C | ±0.2 °C | 2 | (30 x 15.5 x 15 cm) | (48 x 22 x 34 cm) | | | 220–240, 50/60 | 78905-21 |
| 11 L Ambient +5 to 100 °C | Ambient +5 | ±0.2 °C | 4 | 11 ¹³ / ₁₆ x 97/ ₁₆ x 77/ ₈ in (30 x 24 x 20 cm) | 18 ¹⁵ / ₁₆ x 11 ¹³ / ₁₆ x 18 ¹⁵ / ₁₆ in (48 x 30 x 48 cm) | 29 lb (13 kg) | 800 W | 120, 60 | 78905-22 |
| | to 100 °C | | | | | | | 220–240, 50/60 | 78905-23 |
| 27 L | Ambient +5 to 100 °C | ±0.2 °C | 8 | 19 ¹¹ / ₁₆ x 11 ⁵ / ₈ x 7 ⁷ / ₈ in (50 x 29.5 x 20 cm) | 26 ¹³ / ₁₆ x 14 ³ / ₁₆ x 15 ³ / ₈ in (68 x 36 x 39 cm) | 36 lb (16 kg) | 1000 W | 120, 60 | 78905-24 |
| | | | | | | | | 220–240, 50/60 | 78905-25 |

^{*}The 6.6 L models do not accept the 78905-29 tube rack.


Cole-Parmer® WB-200 Series Water Baths

Optional Tube Racks

| Accommodates | Module width | Fits bath size | Item number |
|---------------------------|--------------|-------------------|-------------|
| 20 x 13 mm diameter tubes | 1 | 6.6, 11, and 27 L | 78905-26 |
| 20 x 18 mm diameter tubes | 1 | 6.6, 11, and 27 L | 78905-27 |
| 5 x 31 mm diameter tubes | 1 | 6.6, 11, and 27 L | 78905-28 |
| 8 x 56 mm diameter tubes | 3 | 11 and 27 L | 78905-29 |

Optional Holed Flask Racks

| Accommodates | Fits bath size | Item number |
|---------------------------------|----------------|-------------|
| 2 flasks (33 to 92 mm diameter) | 6.6 L | 78905-30 |
| 4 flasks (33 to 92 mm diameter) | 11 L | 78905-31 |
| 8 flasks (33 to 92 mm diameter) | 27 L | 78905-32 |

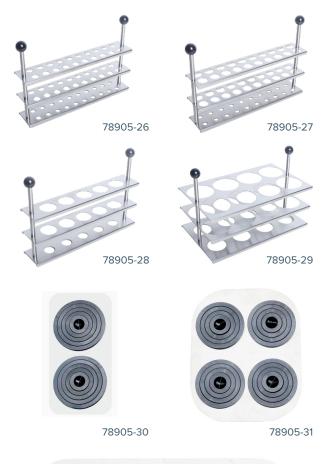



78905-32

coleparmer.com

USA: +1.800.323.4340 +1.847.549.7600 Canada: +1.800.363.5900 China: +86.21.5109.9909 France: +33 (0) 1486 37800 Germany: +49 (0) 9377 92030 India: +1.800.266.1244 Italy: +39 (0) 1313 89513 UK: +44 (0) 1480.272279

All other countries: +1.847.549.7600

