

TFAP4 Polyclonal Antibody

Catalog number: 12017-1-AP


Size: 41 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

TFAP4; AP-4; bHLHc41

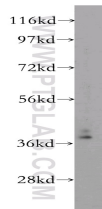

human brain tissue were subjected to SDS PAGE followed by western blot with 12017-1-AP(TFAP4 Antibody) at dilution of 1:500

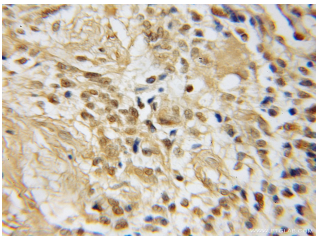

Immunohistochemistry of paraffin-embedded human gliomas using 12017-1-AP(TFAP4 Antibody) at Dilution 1:50 (under 10x lens)

Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	Human brain tissue
Calculated TFAP4 MW:	338aa,39kd
Observed TFAP4 MW:	39kd
Positive IHC detected in:	Human gliomas
Recommended dilution:	WB: 1:200-1:2000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag2633
GenBank accession number:	BC010576
Gene ID (NCBI):	7023
Full name:	Transcription factor AP-4 (activating enhancer binding protein 4)

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.