

Key Features:

- VersaCap® body with advanced ergonomic features for comfortable handling
- VersaCap® technology enables flexible, interchangeable connector options
- Transparency with approximately 30% haze
- Leak-proof cap design
- Round Carboy ideal for mixing
- Made in Class 7 Cleanroom environment

Applications:


- Buffer and media preparation
- Storage and transfer
- Cell culture media
- Sampling and solution mixing

The Foxx Life Sciences EZclear® Polypropylene (PP) Round Carboy is a translucent container that allows easy visual inspection of contents while maintaining the durability, chemical resistance, and autoclavability of standard PP carboys. It delivers the clarity user need at most cost-effective price point, combining visibility with the proven performance of polypropylene.

These Round Carboys are perfect for applications that requires homogeneous mixing. The shape eliminates any pockets of undissolved material that square carboys might develop due to uneven corners. The Carboy features include large top handles, easy-to-read metric graduation marks certified to $\pm 5\%$ accuracy, and the carboy material molded into the side for easy identification.

Carboys are available in 10L, 20L and 50L volume and is manufactured from USP Class VI, polypropylene resin, offering good chemical resistance and thermal stability. The carboys are autoclavable, with a recommended sterilization cycle of 121 °C at 15 psi for 20 minutes.

- 1 Large top handles for a full hand grip
- 2 VersaCap® technology with versatile adaptors and replaceable gasket
- 3 EZfil™ wide-mouth shape for easy filling and cleaning
- 4 Carboy base with 25% curve in the bottom
- 5 Laser etched with clear graduation marks up to $\pm 5\%$ accuracy

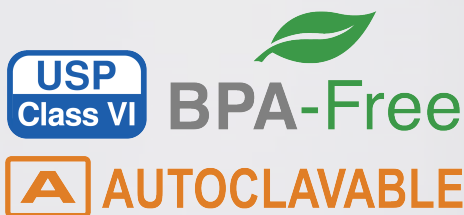

Ordering Information:

Product Number	Description	Volume	Material	Case Qty
515-EL51-FLS	EZclear® Round Carboy 50L, PP, Class VI, 83B Closed VersaCap®, Non-Sterile	50L	PP	1/EA
515-EL52-FLS	EZclear® Round Carboy 50L, PP, Class VI, 83B Closed VersaCap®, Sterile	50L	PP	1/EA
515-EH51-FLS	EZclear® Round Carboy 20L, PP, Class VI, 83B Closed VersaCap®, Non-Sterile	20L	PP	1/EA
515-EH52-FLS	EZclear® Round Carboy 20L, PP, Class VI, 83B Closed VersaCap®, Sterile	20L	PP	1/EA
515-EG51-FLS	EZclear® Round Carboy 10L, PP, Class VI, 83B Closed VersaCap®, Non-Sterile	10L	PP	1/EA
515-EG52-FLS	EZclear® Round Carboy 10L, PP, Class VI, 83B Closed VersaCap®, Sterile	10L	PP	1/EA

Explore our Carboys:

- PP
- HDPE
- PC
- PETG



EZBio® Titanium Round




EZLabpure™ Round




EZgrip® Rectangle


Discover More Single-Use Options:

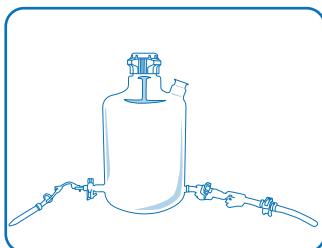

Custom Set

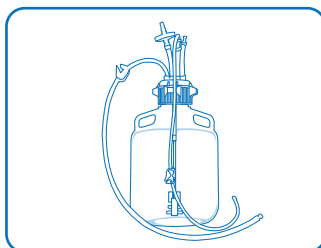

Mixing Set

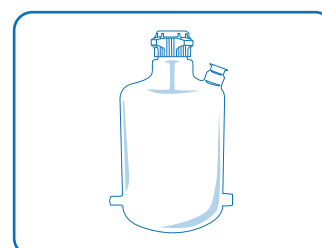

Welded Set

Contact Foxx life Sciences Sales team for more information at sales@foxxlifesciences.com